***Speed* Mam Color™**

English	1
Français	7
Español	13
Português	19
Italiano	25
Ελληνικά	31
Deutsch	37
Netherlands	43
Interpretation	49

SpeedMam Color™

Veterinary diagnostic kit

For *in vitro* use only

■ CLINICAL APPLICATION

Speed Mam Color is a diagnostic tool for the identification of pathogenic bacteria responsible for bovine mastitis and the rapid determination of antibiotic sensitivity. This in-practice test is easy to perform directly on a milk sample. The results of sensitivity profile appear within 24 hours and the identification of the causative germs within 48 hours.

Speed Mam Color is a solution for the individual and specific in-field diagnosis of mastitis: it enables the selection of the most appropriate antibiotic treatment from the vast array currently available to the veterinarian. It takes into account the synergistic or antagonistic effects of the various pathogenic agents, when their concentration is greater than or equal to 10^3 CFU/mL, directly on a sample of mastitic milk.

■ PRINCIPLE

Speed Mam Color is a bacteriologic diagnostic test for farm animals, composed of a culture gallery with:

- 14 wells enabling the determination of the sensitivity profile of the pathogens present in the sample for 14 different antibiotics (or combinations of antibiotics).
- 8 wells for the identification of pathogenic bacteria.
- 2 control wells:
 - the ⊕ bacterial growth control well: a colour change in this well indicates the presence of bacterial concentrations greater than or equal to 10^3 CFU/mL.
 - the ⊖ negative control well: a colour change in this well during the reading time of the test indicates that the test is invalid.

Speed Mam Color is performed from a milk sample. After inoculating the gallery, incubation at a controlled temperature ($+37^{\circ}\text{C}$) enables on the one hand the evaluation of the development of pathogenic microorganisms in the presence of antibiotics commonly used in veterinary practice, and, on the other hand, the identification of the genus of the microorganism(s) present in the sample, via a simple colour change in each well.

■ PROCEDURE

►FOR EACH TEST YOU WILL NEED:

One gallery, one bottle of *Growth medium*, the bottle of *Staph Supplement*, the bottle of *Paraffin oil*, one pipette, one gallery holder, one incubator and one results sheet.

Never use reagents from different boxes.

• PRECAUTIONS:

- It is advisable to wear protective gloves, protective clothing when taking the sample

and during all stages of the test procedure (inoculation and reading).

- In case of contact of one reagent with skin, remove/take off immediately all contaminated clothing and rinse cautiously with water for several minutes.
- In case of contact of one reagent with eyes, rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.
- Perform the test procedures on a dry, clean, and disinfected work surface.

• **SAMPLING:**

- The milk must be sampled in an aseptic manner and stored immediately in the fridge.
- We recommend storing a part of the sample between +2°C and +8°C for 48 hours (the time needed to obtain a complete result) in case further testing is required to confirm results.
- Before proceeding with the test, we recommend allowing the milk sample to come to room temperature for 30 minutes (or warm it gently in an incubator).

• **PREPARING THE GALLERY:**

- Using a gallery taken from its pouch just prior to use, record the name of the animal and the date of culture on the label. Remove the adhesive label covering the gallery.
- Stick the long edge of the label along the long edge of the gallery, to enable access to all of the wells whilst keeping their identification visible.
- Please note that prior to bacteria growth the content of the control wells and of the antibiotic wells may not be visible with the naked eye.

• **PREPARING THE SAMPLE:**

- After homogenization of the milk sample, **add 3 drops of milk to the bottle of *Growth medium*, using the supplied pipette.** Then, homogenize by slowly inverting the bottle.

• **INOCULATING THE GALLERY:**

- **Inoculation** of the wells:

Through the included dropper cap, **add 3 drops of the inoculated *Growth medium* to each well.**

- Add **2 drops of *Staph Supplement* to the STAPH well.**

- Add **2 drops of *Paraffin oil* to each well except the E. COLI, PSEUDO and STAPH wells.**

- Reposition the adhesive label on the gallery and make sure that you remove the transparent film from the back first for improved adherence.

- Place the gallery on the rigid card holder to improve the reading contrast and maintain a homogeneous temperature in all wells during culture.

• **CULTURE:**

- After inoculation, immediately incubate the gallery at +37°C in an incubator.

• **READING:**

1 - Reading the control wells:

Read control wells after 24 hours of incubation at +37°C. (plus/minus 2 hours)

- The \ominus negative control well should remain colourless:

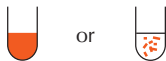
- A change in the \oplus growth control well from colourless to red or with red streaks, indicates that the concentration of pathogenic bacteria is greater than 10^3 CFU/mL

Absence of pathogenic bacteria
(except *Mycoplasma spp.*)
The \oplus well remains colourless

Presence of pathogenic bacteria
The \oplus well turns red or contains red streaks

In this case, the antibiotic sensitivity profile can be read.

- The sample can be considered not infected by pathogenic bacteria (except *Mycoplasma*) if the \oplus well has not changed its colour after 48 hours.

Particular case of *Mycoplasma infections*:

The \oplus well and the antibiotics wells of Speed Mam Color do not allow a growth of *Mycoplasma spp.* The growth of *Mycoplasma spp.* **requires 7 days of incubation at +37°C with Speed Mam Color and can be achieved only into the MYCOP well.**

For these reasons, the \oplus well cannot be an indicator for the presence or the absence of *Mycoplasma spp.* into the sample. **Moreover, *Mycoplasma spp.* does not show any antibiotic susceptibility profile on the gallery.**

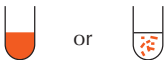
2 - Reading the antibiotic wells:

Read antibiotic wells immediately after interpretation of the control wells

Antibiotic wells can be interpreted only if negative control well is colourless and red color is observed into positive control well.

No colour change.
No bacterial growth.
Bacteria SENSITIVE to the antibiotic

Turns red or presence of red streaks.
Bacterial growth
Bacteria RESISTANT to the antibiotic

Antibiotic table of Speed Mam Color:

The identification of the antibiotic present in each well is printed on the sticky label, above each well.








Abbreviation	Antibiotic	Abbreviation	Antibiotic
CLO	Cloxacillin	GEN	Gentamicin
AMC	Amoxicillin + Clavulanic ac d	SPI	Spiramycin
AMP COL	Amp cillin + Colistin	TYL	Tylosin
CFL	Cefalexin	MAR	Marbofloxacin
CFP	Cefoperazone	PEN DHS	Penicillin + Streptomycin
CFT	Ceftiofur	SUL TMP	Sulfonamides + Trimethoprim
CFQ	Cefquinome	TET NEO BAC	Tetracycline + Neomycin + Bacitracin

3- Interpreting the bacterial identification:

Read identification wells after 48 hours of incubation at +37°C, plus/minus 2 hours – (24h after control and antibiotic wells interpretation)

Only Mycoplasma well can be read after 7 days of incubation at +37°C.

The test may reveal associations of different bacteria.

	Negative control	Positive control	ATB sensitive	ATB resistant	STAPH	STREP	STREP + ESC	ENTERO COCCUS	PSEUDO	E. COLI	ENTERO BACT	MYCOP
<i>Staph. spp.</i>	 Colourless	 Red	 Colourless	 Red	 Yellow	 Colourless	 Colourless	 Colourless	 Colourless	 Colourless	 Red	 Orange
<i>Strep. agalactiae</i>	 Colourless	 Red	 Colourless	 Red	 Red	 Grey-blue	 Colourless	 Colourless	 Colourless	 Colourless	 Red	 Orange
<i>Strep. dysgalactiae</i>	 Colourless	 Red	 Colourless	 Red	 Red	 Grey-blue	 Colourless	 Colourless	 Colourless	 Colourless	 Red	 Orange
<i>Strep. uberis</i>	 Colourless	 Red	 Colourless	 Red	 Red	 Grey-blue	 Black	 Colourless	 Colourless	 Colourless	 Red	 Orange
<i>Enterococcus</i>	 Colourless	 Red	 Colourless	 Red	 Red	 Grey-blue	 Black	 Brown-Black	 Colourless	 Colourless	 Red	 Orange
<i>Pseudomonas</i>	 Colourless	 Red	 Colourless	 Red	 Red	 Colourless	 Colourless	 Colourless	 Navy-blue	 Colourless	 Red	 Orange
<i>E.coli</i>	 Colourless	 Red	 Colourless	 Red	 Red	 Colourless	 Colourless	 Colourless	 Colourless	 Blue-ring	 Orange-Yellow	 Orange
<i>Enterobacteriaceae (Except E.coli)</i>	 Colourless	 Red	 Colourless	 Red	 Red	 Colourless	 Colourless	 Colourless	 Colourless	 Colourless	 Orange-Yellow	 Orange
<i>Mycoplasma</i>	 Colourless	 Colourless	 Colourless	 Colourless	 Red	 Colourless	 Colourless	 Colourless	 Colourless	 Colourless	 Red	 Pink Red

Comments:

Incomplete colour change:

If an identification well turns to an intermediate colour (often associated with turbid medium) after 48 hours of incubation at +37°C (e.g.: STAPH well turns from red to turbid orange), the well should be re-examined again 24h later (at 72 hours after inoculation)

Result at 72 hours:

- Complete colour change: bacterial identification confirmed.
- Intermediate colour persists: bacterial identification is considered negative.

■ INTERPRETATION

Result sheets are supplied to record the identified bacteria and their antibiotic sensitivity and resistance profiles. The choice of antibiotic should take into account the results of the analysis and the pharmacokinetics of the chosen molecule.

■ RECOMMENDATIONS

• STABILITY / STORAGE

- The kit is stable at between +2°C and +8°C for 16 months from the date of manufacture (see expiry date on the label of the kit). Do not expose the kit to temperatures below freezing (0°C).
- It is advisable to leave all of the reagents and the gallery at room temperature for at least 15 minutes prior to use.
- Never use bottles of reagents from different boxes.

• PRECAUTIONS OF USE:

- The volume and the number of drops of sample and *Growth medium* recommended in the protocol are determining for the quality of the results.
- During reading, we recommend leaving the gallery on the rigid card holder.
- If the ⊕ well does not change colour after 24 hours of incubation at +37°C, leave it to incubate for a further 24 hours to confirm the absence of pathogen in the sample (except *Mycoplasma spp.*).

The above recommendations are only guidelines; no diagnostic method is 100% accurate. The purpose of this test is to identify the bacterial agent responsible for bovine mastitis and to assist the veterinarian in his choice of treatment. The interpretation of the test by the veterinarian should always take into account the history, clinical examination, and the results of any further diagnostic tests. The definitive diagnosis is the prerogative and responsibility of the veterinarian.

Bio Veto Test and its distributors cannot be held responsible for any consequences linked to the incorrect usage of this test or misinterpretation of the results.

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